

GENIUS LIGHT

Multiprogram Hi-tech Care

QUICK START GUIDE



The 1st multifunction mask combining Light Therapy® & Electrostimulation.




Beauty on demand

8 beauty programmes, 2 options:

Light mode: LED combined with ionophoresis microcurrent.

Light mode + Electrostimulation: LED combined with different types of electrical stimulation.

A. DEVICE OVERVIEW

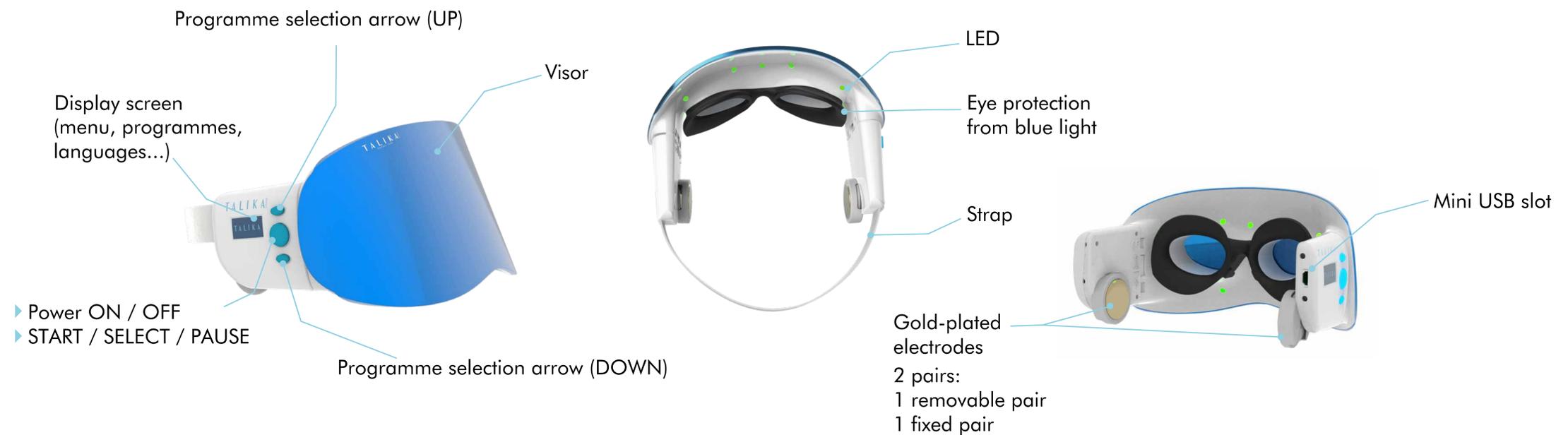


B. BEFORE USE

1. Cleanse your face thoroughly. Do not use **GENIUS LIGHT** while wearing make-up.
2. Tie back your hair.



3. Before switching on **GENIUS LIGHT**, first put it on and adjust the strap length and gold-plated electrodes for optimal and comfortable fit.

If the fit of the electrodes is too tight, take off the removable pair of electrodes and adjust the position of the main electrodes by sliding it. To adjust electrodes position slide them back or forward.

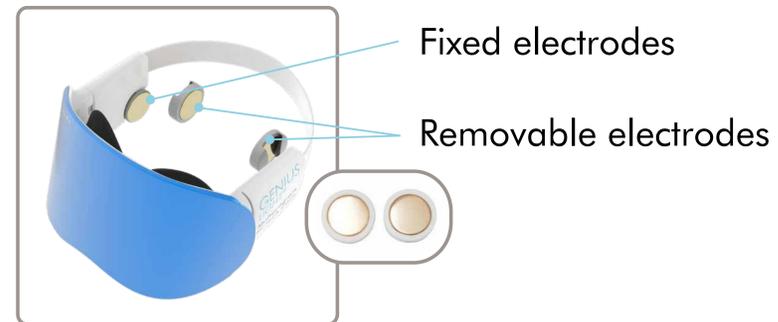

4. Humidify your face using cosmetic sponge or water spray.



Apply water on the area of skin contact with electrodes.

5. Set up **GENIUS LIGHT** beauty care programme.

C. PROGRAMME SET UP

GENIUS LIGHT CONTROL BUTTONS




1. SWITCH ON GENIUS LIGHT

Press START button to turn on GENIUS LIGHT. Display screen will light up with TALIKA GENIUS LIGHT logo.
Press START button once more to go to SELECT LANGUAGE screen.



- UP
- START / SELECT / PAUSE
- DOWN

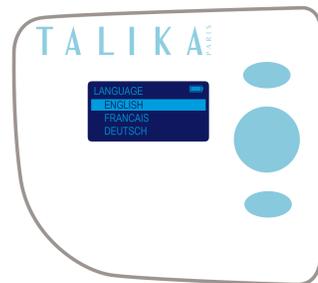

2. SELECT LANGUAGE

Use UP and DOWN buttons to select one of the following languages:

- ENGLISH
- FRANCAIS
- DEUTSCH
- ITALIANO
- ESPAÑOL
- PORTUGES
- NEDERLANDS
- РУССКИЙ
- 中文
- 日本語



- UP
- START / SELECT / PAUSE
- DOWN



3. SELECT MODE

Use UP and DOWN buttons to select MODE:


- ▶ LIGHT
- ▶ LIGHT + Electrostimulation

Press SELECT to confirm.

To return to previous MENU select BACK arrow sign using UP button and press SELECT.



- UP
- START / SELECT / PAUSE
- DOWN
- BACK

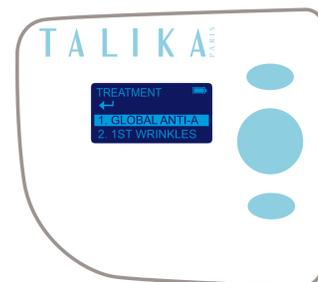

4. SELECT PROGRAMME

Use UP and DOWN buttons to select one of the following programmes:

1. GLOBAL ANTI-AGING
2. 1ST WRINKLES
3. ANTI DARK SPOTS
4. ANTI-BLEMISH
5. DETOX
6. STRESS RELIEF
7. RADIANCE BOOST
8. AFTER SUN

Press SELECT to confirm.

To return to previous MENU select BACK arrow sign using UP button and press SELECT.



- UP
- START / SELECT / PAUSE
- DOWN
- BACK

D. PROGRAMMES



LIGHT MODE



5. START PROGRAMME

- ▶ START icon will appear on the screen after confirmation of your selection. You can begin the care.
- ▶ Make sure your face is still damped. If your face gets dry during the device set up, humidify your face again.
- ▶ Put **GENIUS LIGHT** on before pressing START to start the beauty care in order to avoid being dazzled by the lights.



Ionophoresis is an imperceptible & fixed intensity current. No need to adjust the intensity level. This current does not involve any sensation except from a slight tickle for some types of skin.



The duration of programme is pre-set. The count down clock on the display shows time remaining programme duration.

Programme will stop automatically after pre-set time is run off.



If you don't use the device for more than 3 minutes it will SWITCH OFF automatically.


To force the device to shut down press START button for 2 seconds. The screen will display "BYE" and the device will be shut down.

LIGHT MODE + ELECTROSTIMULATION



5. START PROGRAMME

- ▶ "Ready to Start" message will appear on the screen after confirmation of your selection. You can begin the beauty programme.
- ▶ Make sure your face is still damped. If your face gets dry during the device set up, humidify your face again.
- ▶ Put on **GENIUS LIGHT**. The programme will start automatically with the contact of the gold-plated electrodes on your skin.

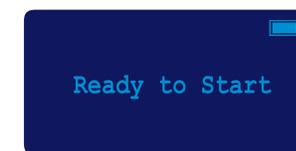

Electrostimulation programme is an electrical stimulation that will make your face muscles move. You can adjust the intensity from 1 to 20, depending on your skin feeling.



There are 20 level of intensity shown on the display next to the count down clock for the programme. Use UP and DOWN buttons to adjust intensity.



IF YOU FEEL ANY UNCOMFORTABLE ELECTRICAL SENSATION: you can either decrease the intensity by scrolling down or press START / SELECT / PAUSE to stop the device.

The duration of programme is pre-set.

Programme will stop automatically after pre-set time is run off.



If you take off the device before the end of programme, it will get on "PAUSE" mode automatically.

If you don't use the device during 3 minutes it will SWITCH OFF automatically.



If you don't use the device for more than 3 minutes it will SWITCH OFF automatically.

To force the device to shut down press START button for 2 seconds. The screen will display "BYE" and the device will be shut down.

E. CHARGING



Take the cable and connect the mini USB part on the mask to a USB charger with adapter or straight into a computer port.



Battery power status is indicated on display screen.



BATTERY FULL

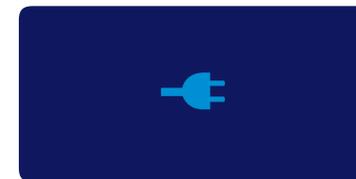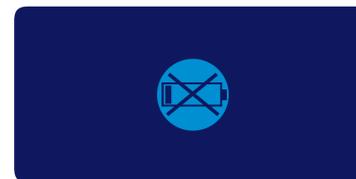

BATTERY NEEDS TO BE CHARGED (ICONS WILL ALTERNATE)



CHARGING

GENIUS LIGHT device cannot work while charging.